

Be sure to read this instruction manual before using this product.  
Keep this manual available for reference when needed.

## Instruction for Use

Trade Name: Micro Scissors

**Warning**

Please read these instructions carefully prior to using this product. Product should be used according to these instructions and pay close attention to the safety of patients. Failure to follow manufacturer's recommendations may cause harm or injury to patient.

**For the US market**

Do not reuse the device when it is used on a patient with, or suspected of having Creutzfeldt-Jakob Disease (CJD) or variant CJD (vCJD).

**For markets outside the US**

When the device is used on a patient with, or suspected of having CJD or vCJD, make sure to adhere to the most recent and updated restrictions available in each country and/or region for its reuse.

Be cautious of the possibility of secondary infection. Refer to [www.a-k-i.org](http://www.a-k-i.org) or AAMI Standard ST79 for more information related to cleaning and sterilization.

**Contraindication/Prohibition**


1. Use for intended purpose only  
Use devices for their intended purposes only. Incorrect use could cause this product to break.
2. Prohibition of use of chemicals  
Avoid exposing this product to chemicals. Doing so could damage this product due to corrosion.
3. Prohibition of secondary processing of this product  
Do not apply any secondary processing to this product, for example, do not apply impacts or vibration markings to the surface of this product. Doing so could break this product.
4. Handle with care  
Handle this product with care, as it can be deformed or damaged. Rough handling could significantly reduce the service life of tools and appliances.
5. Prohibition of use of polishing powder and wire wool  
When cleaning this product, do not attempt to polish its surfaces with rough polishing powder or wire wool. This could cause scratches on the surfaces of this product and result in rust or corrosion.
6. Prohibition of use of household detergents  
Use only medical detergents to clean this product. Do not use any household detergent. Washing this product with an improper detergent could result in discoloration or corrosion.

**Symbol mark for labeling**

MD : Medical Device

**Specifications**

Shape

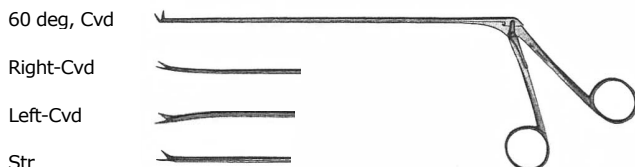

Code No.	Product Description
07-134-01	Micro Dural Scissors A, 60 deg, Cvd
07-134-02	Micro Dural Scissors B, Right, Cvd
07-134-03	Micro Dural Scissors C, Left, Cvd
07-134-04	Micro Dural Scissors D, Str

Material: stainless steel

**Intended purpose**

An instrument, typically used in surgery to cut, e.g. tissue, cloth or sutures.

**Intended user**

This product is to be used by health care professionals, including but not limited to surgeons, nurses and biomedical technicians.

**Instructions for use**

Before using this product, inspect, wash, and sterilize in accordance with these instructions.

**Warning/Caution**

1. Important fundamental cautions
  - Prior to use the device must be sterilized under the standard sterilization conditions or the validated sterilization conditions for which validity has been proven by medical organizations in each country or region.
2. Defect/Adverse event
  - Defect
    - Corrosion or pitting caused by use of chemicals
    - Damage or breakage caused by the corrosion or pitting
  - Adverse event
    - Broken pieces or damaged pieces left in the patient's body.

**Storage/Life**

1. Please store the device in normal ambient temperature areas. Do not store in areas of high humidity where the temperature may dramatically vary causing condensation. Do not store on or near chemicals as the chemicals may cause damage to the device.
2. Service life of this product: 5 years  
(Subject to following manufacturer's specified maintenance, inspection and proper storage requirements.)

**Maintenance/Inspection**

1. Check prior to each use
  - 1.1 Operational and functional checks  
Conduct daily and pre-operation checks of this product to make sure that it functions properly.
  - 1.2 Inspect the instrument for damage, fractures or cracks.
2. Check after each use
  - 2.1 Immediately wash with clean water  
If exposed to bleach or antiseptic solutions, immediately wash: Wash and rinse with clean water immediately and immerse in neutral enzyme detergent to remove any bleach or antiseptic solution, which may contain chlorine or iodine and can damage the instrument. Manually remove contaminated matter by hand or with an ultrasonic-cleaner.
  - 2.2 Further remove any remaining contamination with a plastic brush.

- 2.3 Select a proper detergent for each decontamination method and maintain appropriate density and handling.
- 2.4 Use a soft towel, a plastic brush, or a water gun for cleaning.
- 2.5 To avoid damage, do not use a metal brush, coarse polishing agents, or apply excessive force when handling the device.
- 2.6 Reverse osmosis water is recommended to wash this product.
- 2.7 Only use reverse osmosis water for the final rinse.

It is recommended to use a washer-disinfector for this device. Thermal Disinfection can be used by following the manufacturer's defined parameters.:

Thermal Disinfection Band: 90-93 °C/194.0-199.4°F, 5-10 minutes (A0 value: 3000-12000) (reference EN ISO 15883-1)

- 2.8 Fully dry this product immediately after washing it.  
Do not leave it wet for a longer time than necessary as residual water may damage the instrument.
- 2.9 Use distilled water or reverse osmosis water  
Use distilled water or reverse osmosis water to wash and sterilize this product. Residual chlorine and organic matters in tap water may cause stainings and/or rust and may damage the instrument.
- 2.10 Use a water-based anticorrosive lubricant  
Lubricating oil is completely removed by washing. Do not use without lubricant oil to sliding part, or galling could occur. After washing this product, apply a water-based anticorrosive lubricant prior to sterilization. Do not use this product without applying anticorrosive lubricant on its sliding part. [Galling could occur.]

**3. Sterilization**

Device must be sterilized by users in accordance with validated sterilization procedures, such as an autoclave, that are regulated by medical organizations in each country or region.

Sterilization of the device may be accomplished by steam. The recommended sterilization parameters are as follows;

ISO/TS17665-2	
Temp.	Minimum exposure time
121°C / 250°F*	15 Min
126°C / 259°F*	10 Min
134°C / 273°F*	3 Min

\*According to the Guideline for Disinfection and Sterilization in Healthcare Facilities ,2008 (Update:May 2019)

**For the US market**

FDA has not approved or cleared medical devices, including sterilizers, for the intended use of reducing the infectivity of TSE agents (i.e., prions).

**For markets outside the US**

When the device is used on a patient with or suspected of having CJD or vCJD, make sure to adhere to the most recent and updated restrictions available in each country and/or region for its reuse.

**Maintenance and check by agents**

For safety use of this instrument, conduct periodic checks by the manufacturer or the agent recognized by the manufacturer. Maintenance and check by other agents could cause the adverse events and the decrease of the performance and the function. To schedule the periodic check, contact your local distributor or the manufacturer.

**Packing**

1 piece per pack

**Warranty**

MIZUHO Corporation will repair defective parts of this product without charge for one year from the date of delivery/installment except for cases of damage caused by a third party' s repair, act of nature, improper use or intentional damage. All other warranty terms and conditions are subject to regulations of MIZUHO Corporation.

**Disposal**

This device must be disposed of in accordance with local regulations. Please contact your local distributor for proper disposal.

**Notice**

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority in which the user and/or patient is established.



MIZUHO Corporation  
3-30-13 Hongo, Bunkyo-ku, Tokyo 113-0033, Japan  
<https://www.mizuho.co.jp>