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INDEX

MIZUHO Corporation

3-30-13 Hongo, Bunkyo-ku, Tokyo 113-0033, Japan

TEL(03)3815-3096 FAX(03)3813-5068

<https://www.mizuho.co.jp>

CK07-939-01WW

2023-07 Ver.11

English.....	1	Italiano.....	65
Deutsche.....	21	Español.....	87
Français.....	43	日本語	109

INSTRUCTIONS FOR THE USE OF SUGITA BOOSTER CLIP

Index

- | | |
|---|---------------------------------|
| 1 . Specifications | 10. Warning |
| 2 . Types | 11. Use of MRI |
| 3 . Spring Type | 12. Caution |
| 4 . Action Mechanism | 13. Defects and Adverse Events |
| 5 . Surface Finish | 14. Storage - Care and Handling |
| 6 . Blade Cross-sectional Shape | 15. Validated Cleaning and |
| 7 . Indication for Use | Sterilization Procedure |
| 8 . Directions for Use | - Sterilization - |
| 9 . Contraindications and
Prohibitions | - Storage - |

1. Specifications

Sugita Booster Clip is used to reinforce the closing force of Sugita Aneurysm Clips Permanent Type (standard and large types) as shown in [Table 1].

“Sugita Booster Clip” herein referred to as “clip”.


This clip should be used only in tandem with a previously applied Sugita Aneurysm Clip. **Never apply directly to the cerebral aneurysm or to any blood vessel.** This clip must be placed as close to the spring of the Sugita Aneurysm Clip as possible and must securely grip the proximal part of the clip blades. [Figure 1]

This clip should not be used with clips other than the standard or large Sugita Aneurysm Clips Permanent Type.

[Table 1]

Applicable Type	Prod. Code
Standard Type	07-940-01~09・64~66・68・69 07-943-01~03
Large Type	07-940-10~18・20~25・47~50 60~63・70~73

SUGITA BOOSTER CLIP

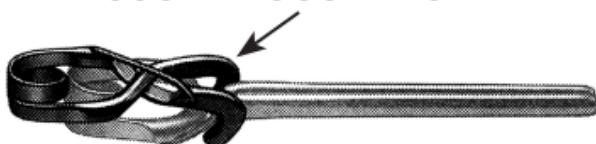[Figure 1]

This clip is made of cobalt-chromium-nickel-molybdenum-iron alloy (ISO5832-7/ASTM F1058 G2).


ASTM: American Society for Testing and Materials

ISO: International Organization for Standardization

This clip is designed for **single-use only**, is exclusively indicated only for use with the Sugita Aneurysm Clips Applying Forceps (hereinafter this is called “clip applier”), and requires cleaning and sterilization prior to use. Reuse may cause breakage of clips and may cause injury or infection.

In order to ensure the traceability, a lot number is shown on each clip.

The closing force of each clip is measured during the final quality assurance stage with a computerized electronic scale at a specified measuring point as shown in **[Figure 2]**. The result of this measurement is shown in newtons and grams on the package label.

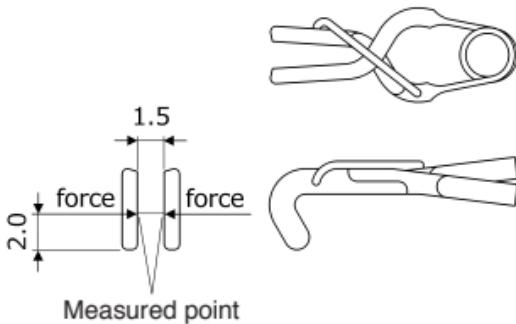[Figure 2]

2. Types

Permanent Type

Type	Color
Booster	Blade tip-yellow

For specific dimensions and shapes, please see this clip's catalog.

3. Spring Type: Integral Coil

4. Action Mechanism: Alligator Action

5. Surface Finish

- Clip Body: Satin
- Material Treatment: Electrolytic Polish

6. Blade Cross-sectional Shape: Rolled

7. Indication for Use

This device is used to reinforce the closing force of standard and large Sugita Aneurysm Clips for permanent occlusion (permanent clips).


8. Directions for Use

This clip is packaged non-sterile, and must be properly sterilized before use.

This clip must be held with the clip applier.

When holding the clip, it must be verified that the spring of the clip is correctly placed in the holding jaws of the clip applier as shown in **[Figure 3]**.

Improper holding of the clip may cause damage or ejection from the clip applier. As this may present a risk during surgery, ensure correct holding of the clip.

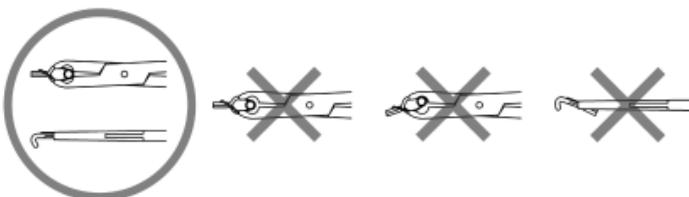[Figure 3]

- Do not hold the clip in angled positions.
 - If the spring of the clip is not properly held with the holding jaws of the clip applier, the clip may slip or spring out from the jaws. This may also damage the
-

spring of the clip, and there is a risk of breakage and reduced holding force.

It is recommended that doctors record the type, lot, and serial number of the clip used in the operation record for each patient.

To support traceability, Patient/Trace labels are included in the clip package.

9. Contraindications and Prohibitions

1) Prohibition of reuse

This clip is for single use only. The clip must be discarded once the blades of the clip have been opened or the clip has come into contact with the blood or tissue of a patient. Reuse may cause breakage of clips and may cause injury or infection.

2) Prohibition of use for unintended purposes

This clip should not be used for purposes other

than reinforcing the closing force of Sugita Aneurysm Clips Permanent Type (standard and large types). Do not use this clip with Fenestrated, Long, or Mini Sugita Aneurysm Clips Permanent Type.

- 3) To prevent changes in the clips closing force or in the mechanical features of the clip when it is handled, only use the exclusive clip applier. Do not use clip applicators from other manufacturers.
 - 4) Do not open the blades of clip directly with fingers. This may cause corrosion or alter the proper holding force. The blades of clip should be opened only when spring of the clip is placed in the proper position inside the clip applier.
 - 5) Modifying or adulteration of this clip is not allowed.
 - 6) In order to avoid a galvanic current, do not use this clip with implants from other manufacturers or with implants made of different materials. A
-

-
- galvanic current may cause corrosion and weaken the holding force.
- 7) This clip must not be used in patients with serious allergies to metals or exogenous materials.
 - 8) This clip must not be used in patients with infectious diseases or sepsis as a complication.
 - 9) This clip must not be sterilized using low-temperature hydrogen peroxide gas plasma.
This clip is not compatible with sterilization using low-temperature hydrogen peroxide gas plasma.
[Sterilization using such method may cause surface color deterioration and affect its characteristics.]

10. Warning

- 1) It is recommended the patient is informed about the procedure and characteristics of this clip.
 - 2) Closing status of the blades of clip should be verified immediately following surgery.
-

[The most serious adverse events that may be reported in conjunction with this clip include displacement of the clip, slip-out of the clip, ejection of the clip from the clip applier, clip fracture, laceration of a cerebral aneurysm, cerebrovascular injury, hemorrhage, death, etc.]

3) Handling of this clip:

Rough handling of this clip may cause deformation or damage to the clip and may harm the patient.

11. Use of MRI

MR Conditional

Non-clinical testing demonstrated that this clip is MR Conditional. A patient with this clip can be scanned

safely in an MR system under the following conditions:

- Static magnetic field of 1.5 tesla or 3 tesla, only
- Maximum spatial gradient magnetic field of 1,000 gauss/cm (10 T/m)
- Maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2 W/kg for 15 minutes of scanning (i.e., per pulse sequence) in the Normal Operating Mode

Under the scan conditions defined, this clip is expected to produce a maximum temperature rise of 2.4 °C (4.32 °F) after 15 minutes of continuous scanning (i.e., per pulse sequence).

In non-clinical testing, the image artifact caused by this clip extends approximately 25 mm from this clip when imaged using a gradient echo pulse sequence and a 3 Tesla MR system.

NOTE: Do not take the clip applier into the MR environment.

The clip applier contains material which has been determined to be MR Unsafe.

12. Caution

Use of this clip is restricted to surgeons with proper training in the use of this clip.

13. Defects and Adverse Events

1) Serious adverse events

- Displacement of the clip (slippage)
- Ejection of the clip from clip applier
- Breakage of the clip
- Rupture of a cerebral aneurysm
- Cerebrovascular injury
- Hemorrhage
- Death

2) Other adverse events

- Infection of surgical wound
- General surgical complications
- Delayed clip displacement due to the application of magnetic field
- Occlusion of adjacent blood vessel
- Foreign body reaction

14. Storage – Care and Handling

This clip is prepared in a non-sterile package and should be stored in its package in a cool, dry, and secure location.

Do not use any clip applier other than the exclusive clip applier.

Handling of this clip should be minimized to reduce contamination or damage.

15. Validated Cleaning and Sterilization Procedure

This clip is supplied non-sterile. Therefore, this clip must be sterilized by the users (healthcare professionals) prior to surgery using validated sterilization processes in accordance with guidelines stipulated by the government.

Proper sterilization and handling of this clip is important for assuring that the clips meet sterile specifications prior to surgery as an aid in minimizing the occurrence of infection in the surgical patient.

Other procedures for minimizing contamination in the operating room, including the use of a surgical unit with laminar flow equipment at the doctor's discretion, are also important.

Sterilization equipment should be validated to assure that it sterilizes specific load configurations properly.

The end user is responsible to validate the sterilization process in the individual sterilizer and assure the sterility.

Sterilization equipment should be cleaned and maintained in accordance with standardized operating procedures specified by its manufacturer.

Clean and sterilize the clip applier according to the clip applier instruction manual.

This clip and instruments must be cleaned to reduce bioburden before first sterilization and resterilization.

Do not use chlorine-based solutions when cleaning this clip. This could cause metallic corrosion.

1) Cleaning

Put the clip in the designated Sugita Clip Case or Sugita Clip Tray (hereinafter this is called “clip case or tray”) and close the lid.

<In case of using Ultrasonic cleaner>

(1) Fully immerse the case (tray) in the water (tap water) of the washing tank.

(2) Ultrasonic cleaning

Cleaning time: 15 minutes or more

Detergent: Neutral enzyme type detergent
(Diluted concentration 1.0%)

(3) Rinsing

- Make sump water with running water, immerse the case (tray) completely in water (tap water or deionized water), shake the case (tray) 5 times up and down.
- Renew the water and shake the case (tray) 5 times up and down in the water again.
- Perform a series of rinsing operations three or more times in total.

(4) Dry thoroughly.

Recommended drying condition: 110 °C (230 °F),
30 minutes or more

<In case of using washer disinfecter (automatic washing)>

Stage	Procedure	Temperature [°C (°F)]	Time [min]	Quality of water	Remarks
I	Preclean	<25 (77)	3	Tap water	—
II	Washing	55 (131)	10	Tap water	Alkaline detergent (such as cleaning agent for aluminum)-Diluted concentration 0.5%
III	Rinsing	<25 (77)	4	Deionized water	Set rinsing process program as "1 time (1 minute) × 4 times = 4 minutes"
IV	Hot water disinfection	—	—	Deionized water	Follow the program of the device.
V	Drying	—	—	—	

2) Inspection

Check clips visually one by one and do not use if the following abnormalities are found.

- Damage
- Deformation
- Corruption
- Inappropriate blade position (deviation, clearance, tip unevenness, etc.)
- Remaining dirt

Since the color of the surface may be faded in repetitive processing, do not use clips that cannot be clearly distinguished with color coding.

3) Packaging

Place clip in the designated clip case or tray and close the lid.

Then pack the case (tray) containing the clip using standard medical packaging material (non-woven fabric or sterilized bag) with protective sterile barrier complying with ISO 11607 for steam sterilization.

- Sterilization -

<Validated Sterilization Method>

Standard followed	ISO 17665-1		
Sterilization method	Pre-vacuum steam sterilization (Autoclave sterilization)		
Sterilization conditions	Sterilization temp.	Retention time	Drying time
	132 °C / 270 °F*	4 minutes or more	20 minutes or more
	134 °C / 273 °F*	3 minutes or more	20 minutes or more

*According to the Guideline for Disinfection and Sterilization in Healthcare Facilities, 2008 (Update: May 2019)

Do not sterilize using low-temperature hydrogen peroxide gas plasma. It may cause discoloration and/or deformation.

With a high-pressure steam sterilizer that has been confirmed to function properly, the validation of sterilization is verified by a medical institution, sterilization conditions with a proven effectiveness shall be applied.

- Storage -

Do not store the clip in high temperatures, or in areas with high humidity and/or where the temperature varies. Please store the clip in clean dust-free and dry area.

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